




SPORICIDAL FOAM

JET is a powerful sporicidal disinfectant that can be used on medical device surfaces located near to the patient.

JET is also a powerful cleaner destroying RNA and DNA in seconds*. Its cleaning efficacy meets the requirements of a medical instrument cleaner, effectively reducing levels of protein and carbohydrate.

JET is the best safeguard for patients and staff wherever the risks of infection are highest. For use in areas such as patient wards, isolation wards, operating theatres, intensive care units, oncology wards and neonatal wards.

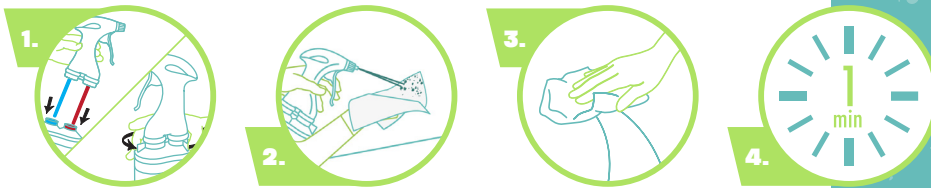

Effective against microorganisms of concern such as:

- *Clostridium difficile*
- Tuberculocidal mycobacteria
- Norovirus
- Coronaviruses**
- Influenza A virus (H1N1)
- Hepatitis B virus (HBV)
- Human immunodeficiency virus (HIV)
- Carbapenem-resistant Enterobacteriaceae (CRE)
- *Klebsiella pneumoniae*
- Methicillin-resistant *Staphylococcus aureus* (MRSA)
- Vancomycin-resistant Enterococci (VRE) *Enterococcus faecium*
- *Candida auris*
- *Escherichia coli*

Complete compliance with EN 14885:2018 and the latest published regulatory standards for disinfectants used in the medical area.

SPORES	EN 17126	YEASTS/FUNGI	EN 13624
	EN 13704		EN 14562
MYCOBACTERIA	EN 14348	BACTERIA	EN 16615
	EN 14563		EN 13727
VIRUS	EN 14476		EN 14561
	EN 16615		EN 16615



**1 UNIFORM CONTACT TIME
OF ONE MINUTE**

**1ST LEVEL OF COMPLIANCE
TO ALL APPLICABLE
EN STANDARDS**

**1 ACTIVE INGREDIENT
CHLORINE DIOXIDE**

**1 DISINFECTANT
FOR ALL MEDICAL
DEVICE SURFACES**

POWERED BY CHLORINE DIOXIDE FOAM

JET is built upon the same core chlorine dioxide technology as the Tristel medical device disinfectant products, but specially adapted for surface disinfection where cleaning and biofilm destruction are so important.

Chlorine dioxide is widely regarded as one of the most effective disinfectants for biofilm removal and prevention.

JET delivers chlorine dioxide as a foam. It creates no aerosols when used.

SAVING SPACE, REDUCING WASTE




3.5 Kilos



0.8 Kilos

Efficacy with the standard wipe in EN 16615 means JET can be used with any sustainably sourced wipe, avoiding the use and waste of plastic wipes.

NEAR THE PATIENT WHERE JET WILL PROVIDE ITS GREATEST PROTECTION:

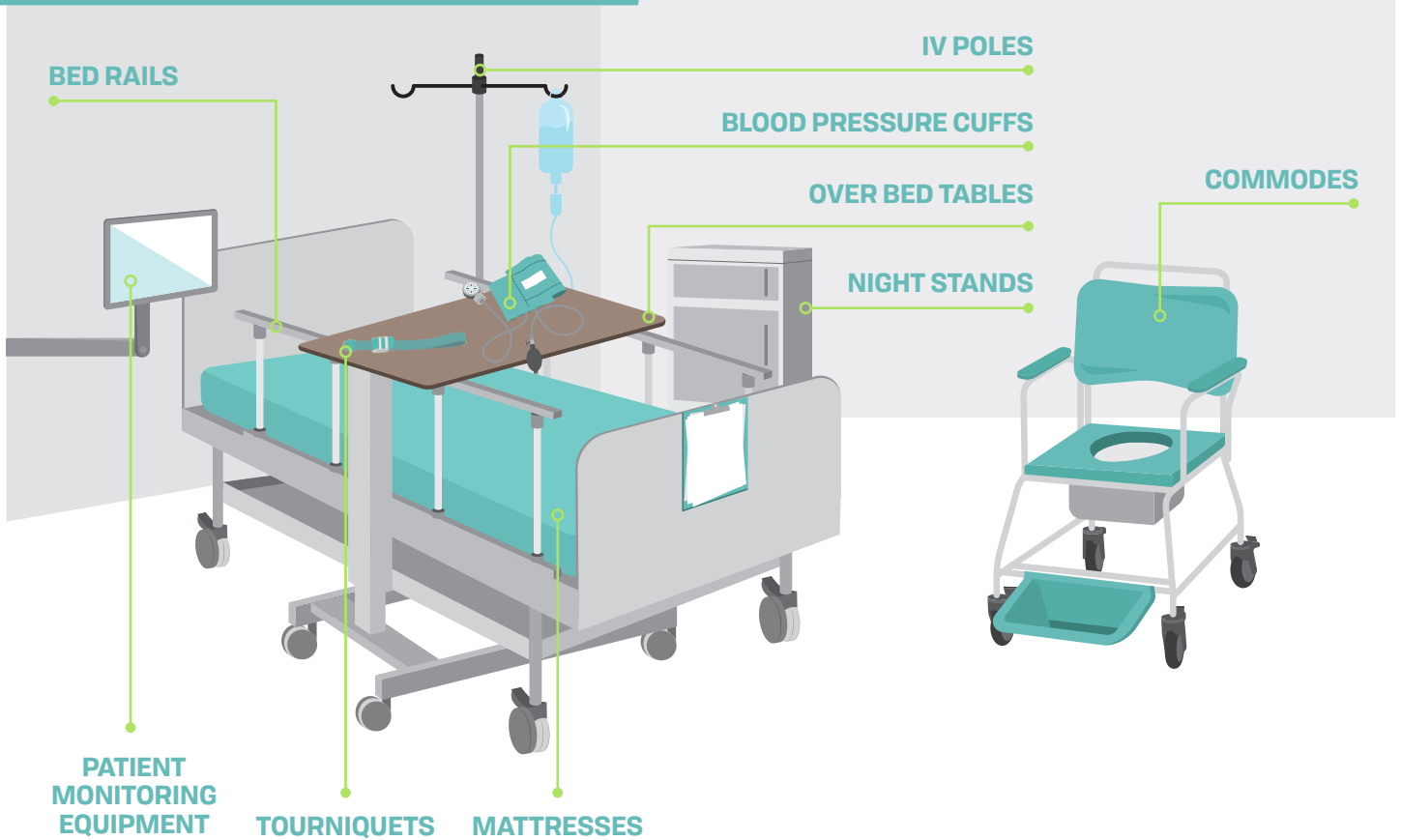

Training and certification is an essential part of customer service with Cache. Ask about our [Online Training Portal](#) for product training at a time and place convenient to you.

**AVAILABLE TO
ORDER NOW!**



2 x 800ml JET per box

Quote: CCH010101

JET is CE marked as a Class IIa Medical Device in accordance with the European Medical Devices Directive 93/42/EEC and the 2007/47/EC amendments thereto.

cache™

Manufactured in the United Kingdom by: Tristel Solutions Limited, Lynx Business Park, Cambs, UK, CB8 7NY
T +44 (0) 1638 721500 - E mail@tristel.com - W www.tristel.com

Hong Kong & Taiwan: Tristel Asia Limited, 21st Floor, 168 Electric Road, Fortress Hill, Hong Kong
T +852 2895 6968 - F +852 2869 4388 - E customerservicehk@tristel.com

Malaysia: Tristel Malaysia Sdn. Bhd, A-25-3A Pinnacle PJ, Lorong Utara C, 46200 Petaling Jaya, Selangor Darul Ehsan, Malaysia
T +603 79317791 - F +603 79319971 - E mail-my@tristel.com

BROUGHT TO YOU BY TRISTEL

Copyright © Tristel Solutions
MkA-80-758-5
October 2020

*Pacific Biolabs (2017) Analysis of Tristel Nucleic Acid degradation properties. **Based on virucidal efficacy validated in accordance with EN 14476 - For Tristel patent information please visit: <http://www.our-patents.info/tristel>